

## INSTALLATION GUIDANCE

# EASY bath seat without backrest

B&M SKU: ID15501.2    Idral code: 15501.2

## Installation sequence

- 1 Assemble the components in the sequence shown in the accompanying Idral assembly drawings.
- 2 Adjust the bath supports within the 500-635 mm range shown in the drawing.
- 3 Before each use, check that the brackets rest correctly at all four contact points on the bath edge. The bath must be stable; take particular care with acrylic baths.
- 4 The bath side walls must have an inclination of at least 5 degrees. Avoid jerky movements while using the seat.

## Important

Maximum load: 110 kg. Check that the bath edges support all four contact points securely before each use. Do not exceed the maximum load.

## INSTALLATION GUIDANCE

# EASY bath seat without backrest


B&M SKU: ID15501.2 Idral code: 15501.2

## Idral assembly details



**PZ2**

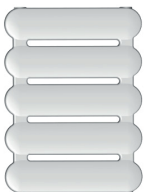


| 1                                                                                   | 2                                                                                   | 3                                                                                   | 4                                                                                   | 5                                                                                    | 6                                                                                     | 7                                                                                     |
|-------------------------------------------------------------------------------------|-------------------------------------------------------------------------------------|-------------------------------------------------------------------------------------|-------------------------------------------------------------------------------------|--------------------------------------------------------------------------------------|---------------------------------------------------------------------------------------|---------------------------------------------------------------------------------------|
|  |  |  |  |  |  |  |
| 1x                                                                                  | 4x                                                                                  | 2x                                                                                  | 2x                                                                                  | 4x                                                                                   | 4x                                                                                    | 12x                                                                                   |

## INSTALLATION GUIDANCE

# EASY bath seat without backrest

B&M SKU: ID15501.2 Idral code: 15501.2

## Idral assembly details

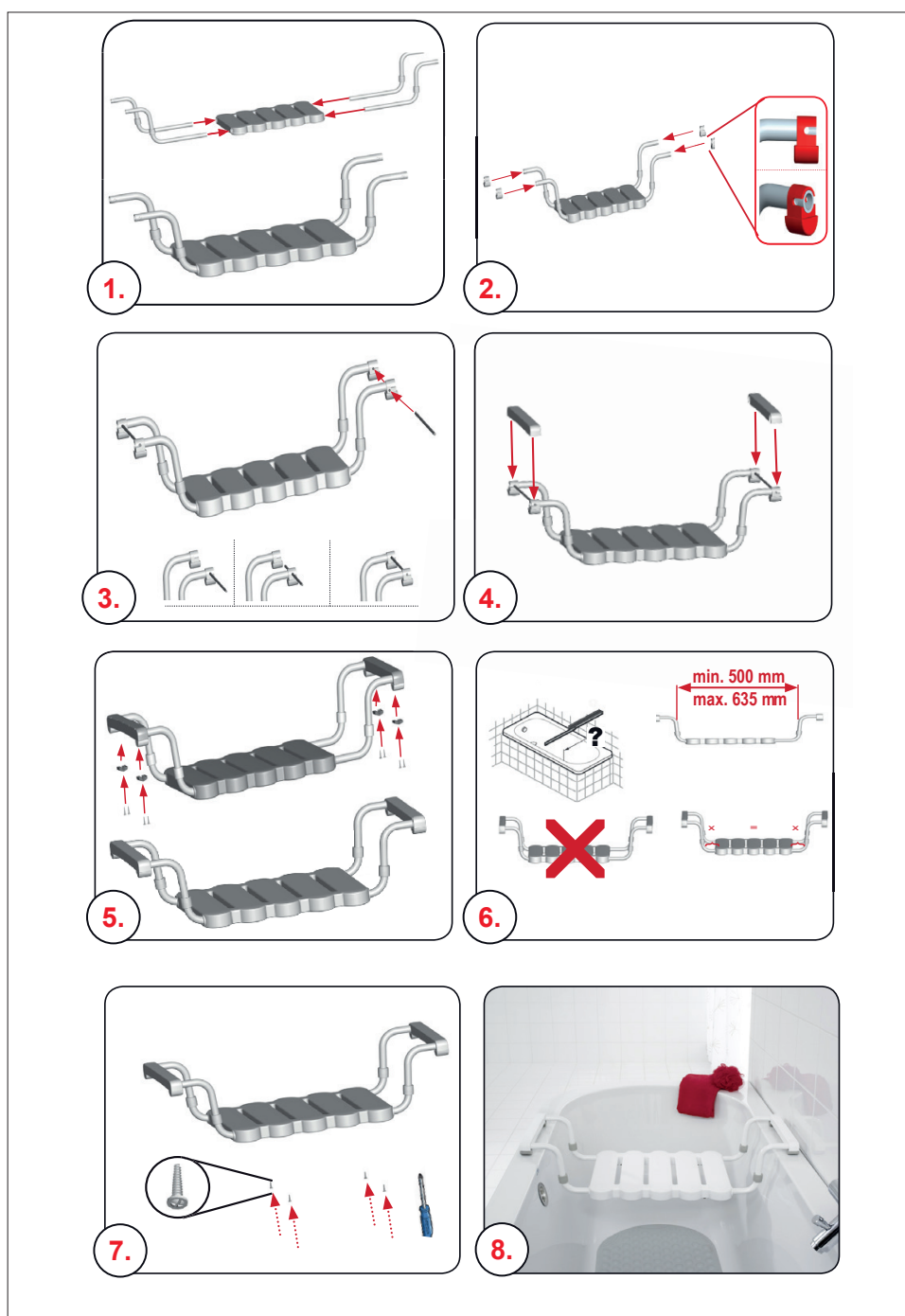

Source: Idral ISTA22-01 / page 7